## D-mine\_FolderPen\_INT\_V03

Non-use starts are started version in the start on contrast on the start on contrast on the start on the star ntho

The characteristic is the characteristic



## D-*mine*® Pen

The new reusable D-mine® Pen injector. This high quality Pen is developed for intermittent bolus injections with Dacepton<sup>®</sup> (Apomorphine 10 mg/ml in 3 ml cartridges) in Parkinson's Disease. The automatic delivery-system enables low and constant injection force for PD patients with motoric impairment. Possible dose correction and last dose stop guarantee safe comfort.

The D-*mine<sup>®</sup>* Pen injection may deliver a controlled "on" period for patients and has therefore a positive impact on daily activities!

## Dacepton



## GET YOUR PATIENTS "ON"

The new D-mine® Pen for subcutaneos intermittent bolus injections


- Reusable Pen injector for Dacepton<sup>®</sup> (Apomorphine 10 mg/ml in 3 ml cartridges)
- Dosing mg-based in 0.5 mg steps
- Automatic delivery-system: constant-spring driven injection for low and constant injection force
- In-use stability for 15 days
- Last dose stop the dialed dose cannot exceed the remaining volume



- Dose correction-module for prevention of faulty doses
- Real time display for interrupted injection
- Audible feedback during dosing, dose correction and delivery
- Intuitive cartridge exchange
- Carrying case for Pen with space for one spare cartridge and needles and comfortable removal of used needles





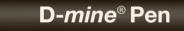



EVER Neuro Pharma GmbH Oberburgau 3, 4866 Unterach/Austria www.everpharma.com, www.d-minecare.com

DAC/INT/06/2021/22